


## Negative Pressure Wound Therapy (NPWT) Dressing Kit

The Cork Medical Products Wound Dressing Kit is designed to be used in conjunction with an NPWT pump to provide negative pressure wound therapy to aid and promote wound healing by removing excess exudates, infectious material and tissue debris.

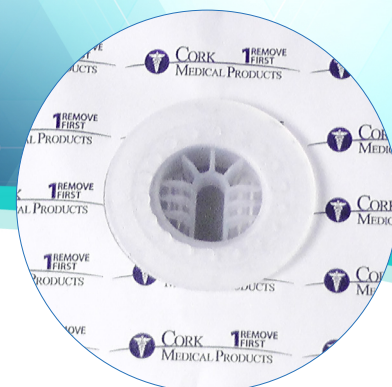

**Advanced Wound Dressing Kit Provides Optional Seal Without Sacrificing Patient Comfort**

- ✓ **1.2.3. film drape with instructions for a quick and user-friendly application**
- ✓ **Soft and comfortable port pad designed to withstand collapsing under NPWT**

### NPWT WOUND DRESSING KIT FEATURES

The Cork NPWT Wound Dressing Kit is comprised of all components necessary to dress a wound for negative pressure wound therapy. The kit's advanced design provides a user-friendly application for patient comfort and allows physicians the best treatment modality for NPWT.

- **1.2.3. film drape** provides clinical breathability, excellent adherence and has step-by-step instructions printed directly on film for a quick and user-friendly application
- **High tensile strength foam** is durable, flexible, and can be cut to contour a variety of wound shapes and sizes, allowing for an individualized fit
- **Dome-like silicone port pad** design is durable enough to withstand collapsing under NPWT yet soft and comfortable
- **NPWT application instructions** are included on the outer pouch to allow for user ease and help eliminate user error

## COMPONENT SPECIFICATIONS

Components within the kit shall be used to dress a wound in conjunction with negative pressure wound therapy per the instructions of a prescribing physician. The kit shall connect to an NPWT pump via the connection tubing provided within the kit.

The components included are as follows:

- Available in three (3) sizes: Reticulated, 30 pores per inch (PPI) polyether, polyurethane foam:
  - Medium foam: 3/4" x 5" x 8" (1); Thick foam (non-insurance coded): 1.25" x 5" x 8" (1); Large foam: 3/4" x 6" x 10" (1)
- Two (2) transparent polyurethane film drapes with adhesive backing, 10" x 12"
- Port pad assembly is comprised of soft silicone port pad, tubing, luer connector and pinch clamp
- One (1) paper ruler to measure wound
- Each kit is single-use and housed in a Tyvek peel pouch; inner kit components are individually pouched and sterilized by Ethylene Oxide or Gamma irradiation
- Medicare HCPCS Code: A6550



In Partnership with

**CORK**  
MEDICAL



For more information, contact the InfuSystem NPWT Team at  
1-844-534-1997 or [npwt@infusystem.com](mailto:npwt@infusystem.com)

